




ICON INSTRUMENTS COMPANY

(Manufacturer of Laboratory Equipment)

IIC-137

WILLY MILL (ARTHUR H. THOMAS TYPE)

- Feature: MS/ CRC Sheet Cast steel chamber duly painted, recommended for milling for sample, animal hair, bones cotton seeds organic drug, Leather, tobacco Wheat & other material
- Designed for preparation of material for lab. Analyses without loss of moisture due to heating.
- With 2 Stationery blade & rotor.
- With 4 cutting edges.
- Sample introduces through hopper & Cover fitted with 1/4HP Motor.
- Three deliveries S. Steel sieves of 20, 40 & 60 mesh with plastic receiver.
- Removable glass plate forms the face of grinding chamber & Permit observation during operation.
- A jar placed at the lower & receives the milled powder.

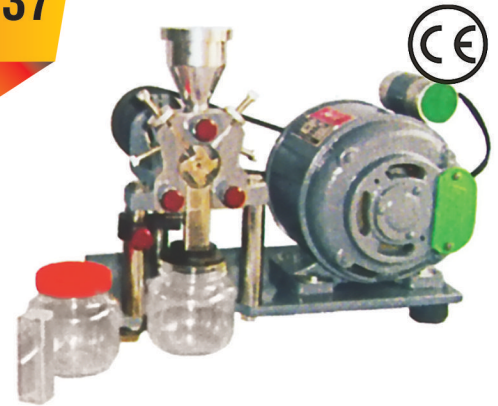

Applications :

For Pharmaceutical Cotton Seeds, Leather & Wheat Industries

Model	IIC 137		IIC 137 A		IIC 137B	
Capacity mm	040 x 25mm		065 x 25mm		100 x 50mm	
Chamber	Cast Steel	S. Steel	Cast Steel	S. Steel	Cast Steel	S. Steel
Motor	1/4 HP		1/4 HP		1/4 HP	

Technical Parameters :

- Material : Cast Steel / S. Steel duly machined
 Power : 220/230V AC supply 50 HZ Single Phase
 Motor Capacity : 1/4HP heavy duty
 Duty Cycle : Continuous

VDRL ROTATOR (CLINICAL)

- Rubber Covered metal platform.
- Platform size 13" x 13" can take cavity slides.
- Platform can take blood bottles & flask by using springs holder / clamps.
- Fixed oscillation 180 RPM per minute.
- Suitable for VDRL Test.

IIC-138


Applications :

For Hospital, Pharmaceuticals, Medical Institute, Research Lab.

Model	IIC 138	IIC 138 A
Type	Fixed Speed 180 RPM	Variable Speed 70-210 RPM
Timer	With Timer 0-60Minute	With Timer 0-60 Minute

Technical Parameters :

- Material : CRC/MS Sheet metal duly painted
 Power Consumption : 220/230V AC Supply Single Phase 50 HZ
 Duty Cycle : Continuous

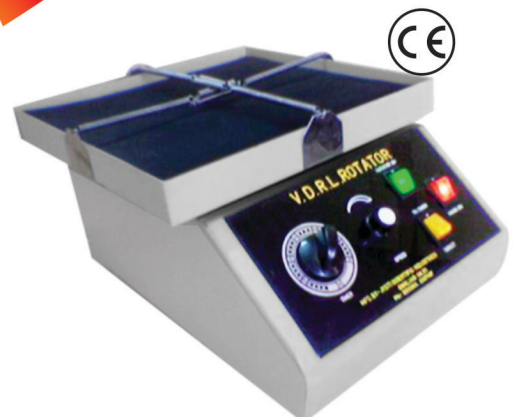

Works : Z-67, Okhla Industrial Area, Phase-II, New Delhi - 110020

Phone : +91-11-26383676 • E-mail : iconinstruments.india@gmail.com, sales@iconinstruments.co.in

Website: www.iconinstruments.co.in